Collection of large amounts of particulate matter from vehicle exhaust using augmented high-volume samplers

Beranek V¹, Stolcpartova J² and Vojtisek M¹

1 Institute for Automobile, Combustion Engine and Railway Engineering, Czech Technical University in Prague 2 Institute of Experimental Medicine, Academy of Sciences of the Czech Republic Contacts -> email: jitka.stolcpartova@biomed.cas.cz

Toxicological aspect of car particle emission is very


Toxicological assays require collection of a rather large amounts of particulate matter: ~0.01 mg for gravimetry ~1 mg for PAH analysis ~100 mg for RNA microarray

Time-consuming and expensive PM sampling

Atmospheric Hi-Vol samplers (2 x 67.8 m3/h) sampling from a full-flow dilution tunnel but this is still not enough & filter capacity only approx. 10 mg due to small, fractal soot particles.

Goal: Innovation of sampling procedure -> collection of higher amount of PM

Two high-volume atmospheric samplers have been augmented with an **auxiliary three-stage blower** (17 kPa @ 2*68 m3/h) and with a **high-strength filter holder** (15 kPa pressure drop ~75 kg).

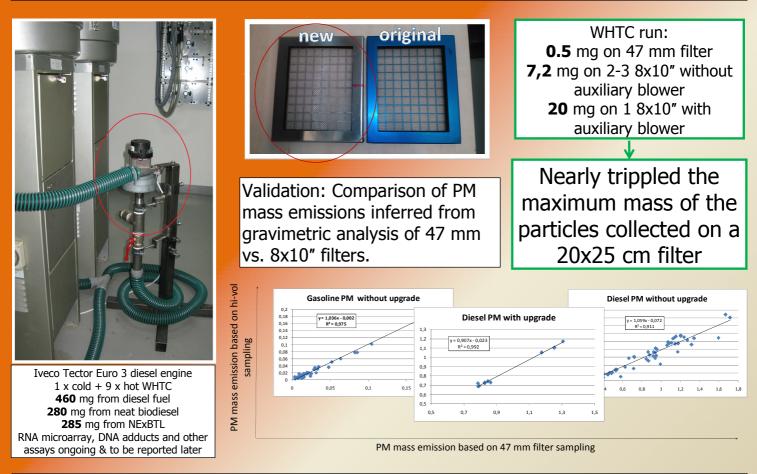

Conclusion

Upgrade of high-volume samplers nearly tripled the maximum mass of the particles collected on a 20x25 cm filter, allowing lower dilution rates, longer tests to be run and more particles to be collected in a given time on smaller number of filters, reducing extraction and analysis cost and time and reducing background effects. Regulation of sample flow at high filter pressure drops still needs to be addressed.

Acknowledgement: The work was funded by the Czech Science Foundation, project BIOTOX - Mechanisms of toxicity of biofuel particulate emissions (13-01438S).